


LIBERTY[®] X

GLIDE
WEAR[®]

Available with GlideWear[®]
Technology. See pg. 25 for details.



*Shown with GlideWear[®]

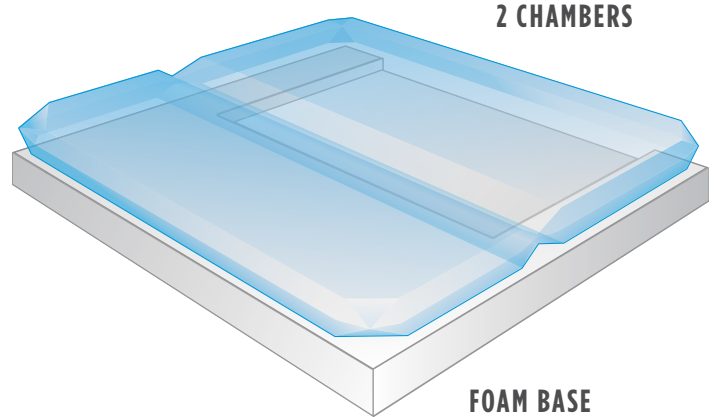

E2622**ADJUSTABLE SKIN PROTECTION****PRODUCT FEATURES**

- The main cushion component is designed with a completely sealed PVC liner containing a front and rear chamber. The cushion liner holds Comfort Cells™ coated in a silicone lubricant giving them a slow moving yet fluid like motion when conforming to the user.
- The liner also has an air valve, that when opened, allows the excess air to escape after a user sits. Once the air has escaped and the cushion has taken on the fitted shape of the user, the air valve can be closed and the cushion will then retain its shape. This vacuum forming technique with the silicone coated Comfort Cells™ captures the user's shape every time.
- A foam base is installed separately under the main cushion component to provide additional femoral support and trochanter loading to the main Vicair® insert.
- A separate cover made of Comfort-Tek™ material covers the foam base giving it superior incontinence protection and infection control. Foam must be removed while laundering cover.

DIMENSIONS

Seat Thickness	4"
Inferior Thigh Support Thickness	4"
Pelvic Contour Height	.5"

STANDARD SIZES

	16" W	17" W	18" W	20" W
16" D	16 x 16		18 x 16	
17" D		17 x 17		
18" D	16 x 18		18 x 18	20 x 18
20" D			18 x 20	20 x 20

Weight Capacity for All Sizes = 330 lbs

LIBERTY® X PART NUMBERS

Cover Style	Stretch-Air™	Comfort-Tek™
Standard	LBX-S-WWDD	LBX-F-WWDD
+ GlideWear®	LBX-SQ-WWDD	LBX-FQ-WWDD

WWDD to be replaced / **WW** = width & **DD** = depth in inches

*Liberty X is not recommended for use in altitudes above 3000 ft.

CLINICIAN TIP

Features a soft feel for those with hypersensitivity

